Power wheelchairs

## Power standing

-5-



### **Corpus VS**

The second second second second

The de march des

Fully programmable standing sequences allow user position to be optimized using supine, semi-reclined, or sit-to-stand techniques. Further customize standing positions with adjustable chest support, multiple knee support options, and power articulating footplates. From face-to-face interactions to improved health outcomes, the Corpus VS seating system empowers you to stand with confidence.

# Designed for independence

New swing-away knee supports are 75% lighter<sup>\*</sup> and feature improved comfort and flexible on-chair storage.

\*As compared to original single-post design

Release mechanism designed for a wide variety of hand function Seat-rail mounted parking spot for convenient onchair storage Integrated BOA® Fit System for a precision fit

P

(r)

BOA is a trademark of BOA.



## Stand without compromise


The F5 Corpus VS combines all the performance features of the F5 Corpus front-wheel drive base with superior power standing functionality. With ComfortRide Pro<sup>™</sup> suspension, experience all-wheel independent, fully adjustable, oil dampened shocks for the ultimate in ride comfort while optimizing traction and performance.

#### F5 Corpus VS Group 4 only:

2-year warranty Preventative Maintenance Plan Hydrographic color options included One-time shroud color change included

Standing & **14**" ActiveHeight

Low, functional seatto-floor heights

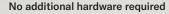



### **Programming simplified**

Programming changes are faster and more intuitive with QuickConfig, the new wireless programming app for Permobil wheelchairs equipped with Power Platform. Designed with clinicians and equipment providers in mind, in-app visuals ease customization of standing sequences, memory positions, and drive performance. The result? Wheelchair setup has never been easier.



Scan to learn more



Available to download from the App store and Google Play

Google Play

Available on the

App Store

#### **Specifications**

	M Corpus VS	F5 Corpus VS
Maximum user weight	265 lb	300 lb
Maximum speed	6 mph, Group 3	5.8 mph, Group 3   7.5 mph, Group 4
Range	20 mi*	20 mi <sup>*</sup>
Base width	25"	25.5"
Base length	36.5"	43"
Base turning radius	22"	30"
Ground clearance	3"	3"
Recommended battery type	Group 24	Group 24
Weight including batteries	424 lb	432 lb
Drive electronics	R-net Power Platform 120A	R-net Power Platform 120A
Seat-to-floor height	17.5", 18.5", 19.5" <sup>†</sup>	17.5", 18.5", 19.5" <sup>†</sup>
Seat elevation	12"	14"
Drive wheel	14"	14"
Casters	7"	8"
Tilt options	Posterior: 0° – 50°	Posterior: 0° – 50°
Power standing	Up to 70°	Up to 80°
Legrest angle - power	85° – 180°	85° – 180°
Backrest angle - power recline	85° – 180°	85° – 180°
Backrest height	20", 23" – 28" (1" incr.)	23" – 28" (1" incr.)
Seat depth	14" – 22" (1" incr.)	14" – 22" (1" incr.)
Seat width	17" – 23" (2" incr.)	17" – 23" (2" incr.)

\* Actual driving range will vary based on driving and battery conditions. |<sup>†</sup> Measured without cushion.

Specifications may vary based upon configuration. Please refer to your owners manual for warnings and instructions for the safe operation of your wheelchair in a variety of driving conditions. Chair weights are approximate and include base, batteries, power seat functions, and seating system without accessories.

App Store is a trademark of Apple Inc. Google Play is a trademark of Google LLC.

